

# **ELECTRO HYDRAULIC OPERATING TABLE**



**Product Catalogue** 



## **Functional Display of Remote Panel**



Table up



Table down



Trendelendburg position



Reverse Trendelendburg position




Turn left



Turn right



Back plate up

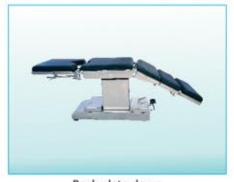

Back plate down



Flex position



Re-flex position



Longitudinal shift towards leg end



Longitudinal shift towards Head end




### STAR-1



#### **General Surgery**

The stability of STAR-1 is suitable for abdominal, chest, and advance laparoscopic surgeries.

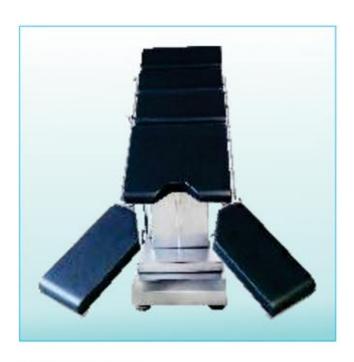

#### **Bariatric Surgery**

Semi standing position achieved by High reverse trendelenberg. Optional width extenders available




Extracorporeal membrane oxygenation (ECMO procedure) with fowlers position, the patient lies on the slightly tilted surgical table. The trendelenburg angle creates respiratory efficiency for patients.



Spine Surgery
Using Laminectomy frame Optional



Lower position for ENT & neurosurgery

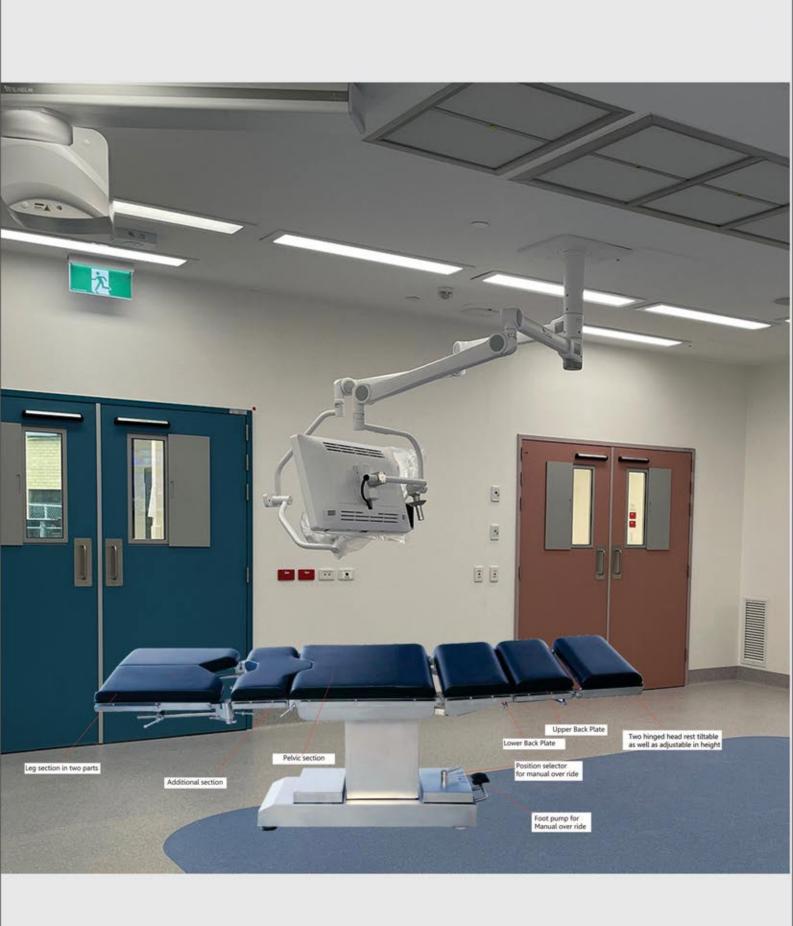



Struma position



Split leg sectionup & down split 180°



Lithotomy position for gynecology & urology



Inbuilt kidney bridge (C-ARM compatible)
There is no bar across the table top


# **Optional Attachment**



Orthotec-2 Radiolucent pelvisplate orthopedic extension



Knee-chest position

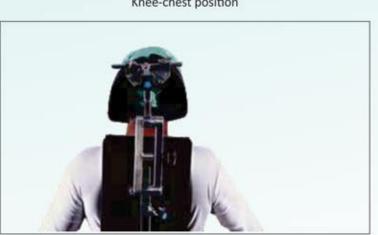

**Shoulder Surgery Attachment** 



Accessories stand



T.H.R. Attachment

### STAR-1

















### **Technical Specifications**

Elevation Range	Highest	1016	From the floor surface to the upper surface of the table top
	Lowest	711°	
Lengthwise Tilt Angle	Head Down	25°	From the level position
	Head up	25°	
Lateral Tilt	Right Down	20°	From the level position
	Left Down	20°	
Back Section	Up	90°	From the level position
	Down	40°	
Leg Plate Flex Angle	Down	100"	From the level position
Sliding table top	,		
Ontional Height A	djustment : 660mr	m - 965mm 7	00mm - 1075mm

#### **Standard Accessories**

- Anesthetic Screen.
- · Shoulder Supports with pads. Arm Boards with pads. Lateral
- Support with pad. Knee Crutches Goepel type. Water proof rubber mattress. Additional Section provided. Body Strap
- X-ray Cassette holder

### **Optional Features**

- The remote is ergonomically designed to have an easy and better grip with spiral cable and have LED backlit screen with symbolic position figures making it convenient to use even in the dark.
- Zero position facility i.e. by pressing the single button OT Table come to normal position.
- Patient reverse orientation mode when head and leg section are interchanged.
- Store up to two preset table top position in its memory which can be recalled anytime by simply pressing M1 or M2 button on the remote.
- Remote control have function of locking operating position to prevent accidental movement of that position during surgery.



### TRINETRA METAL TECH PVT LTD

30 Hari Krishna, Industrial Area, G.T. Karnal Road, Kundli, Sonipat, Haryana-131028

Ph : +91-8802923400

E-mail: trinetrametaltechpvtltd@gmail.com
Web: www.trinetrametaltechpvtltd.com